



...Vero
Senti
 is more than MAGIC!

Play Audiometry

NEW
 Senti Desktop (Vero)



Senti is flexible and easy to operate


Pure Tone Audiometry is still the gold standard for assessing hearing loss. In SENTI, it is available as class 3 audiometer (IEC 60645-1) with additional features: talk over, external input, high frequencies up to 16 kHz...

Speech recognition in noise and speech processing abilities are key skills for children and adults to interact with their environment. SENTI features a speech in noise screening test, image based speech tests for children, screening tests for auditory processing disorders and much more!

The modular approach enables you to assemble yourself the best test battery in your **Senti**

It provides perfect matches for occupational health, general audiology, pediatrics as well as for a doctor's office.

New ideas in audiometry

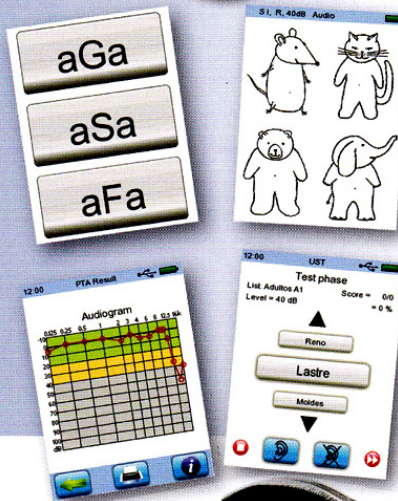

Product Features

- + Large and easy to use touch display of 5.0"
- + Stable desktop housing
- + USB interface to PC
- + External input modality for speech tests e.g for CD
- + Internal microphone for environmental noise recordings, optimisation of measurements
- + Multiple languages support
- + Online help and information available on the device
- + Database to store 1000 patients with multiple tests per patient
- + Direct printing to label printer or via connectivity to PC
- + Updateable any time

Test Modules

- + Pure tone audiometry from 250 Hz up to 16 kHz
- + Levels from -10 to 110 dB HL in 5dB steps
- + Multiple options for accessories: HDA280, HDA200, DD45 (similar to TDH39), B71, free field, insert earphones, patient response switch
- + Bone and air conduction with masking
- + Pure tone audiometry (class 3 EN60645-1)
- + Free-Field option available

Image Based Audiometry

- + Image based audiometry and screening (MAGIC) for children > 4 years

Speech tests:

- + Multiple Speech tests (SUN, MATCH, MAUS, UST...)
- + Screening tests

Distributed by:

HEALTHWORLD INTERNATIONAL PHILIPPINES CORPORATION

0421 Quirino Avenue La Huerta, Parañaque City, Philippines

Telephone: USA +810 4962007 Philippines +632 8290175 Fax: +632 8266770 email: info@healthwordintl.com