



HEALTHWORLD INTERNATIONAL PHILS. CORP.

World Class *QUALITY* Need Not Be Expensive

FDA PHILIPPINES-REGISTERED

...Desktop
Senti
is more than MAGIC!




Senti is flexible and easy to operate

Pure Tone Audiometry is still the gold standard for assessing hearing loss. In Senti, it is available as class 3 audiometer (IEC 60645-1) with additional features: talk over, external input, high frequencies up to 16 kHz...

Speech recognition in noise and speech processing abilities are key skills for children and adults to interact with their environment. Senti features a speech in noise screening test, image based speech tests for children, screening tests for auditory processing disorders and much more!

The modular approach enables you to assemble yourself the best test battery in your **Senti**

It provides perfect matches for occupational health, general audiology, pediatrics as well as for a doctor's office.

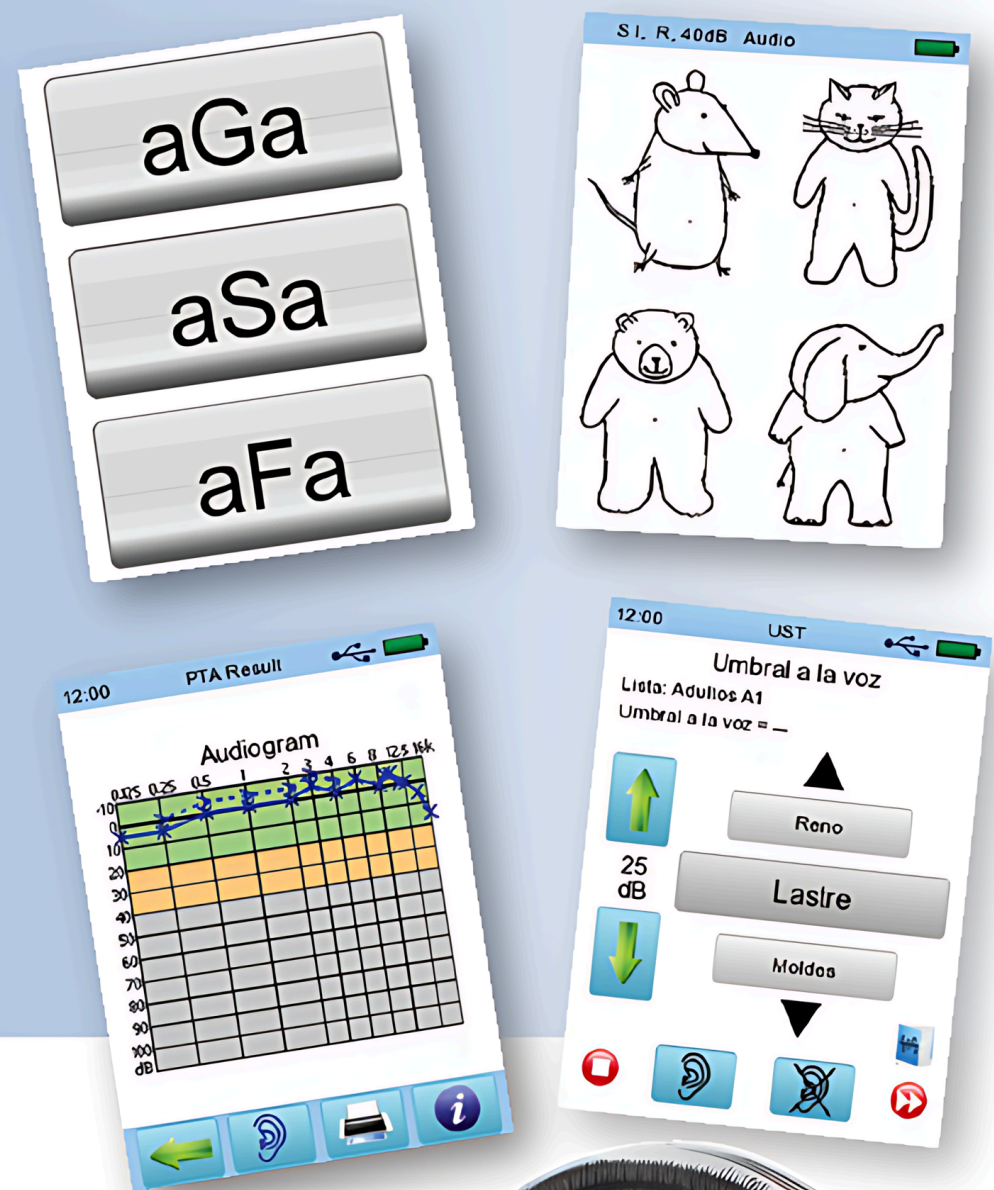

New ideas in audiometry



HEALTHWORLD INTERNATIONAL PHILS. CORP.

World Class *QUALITY* Need Not Be Expensive

...Vero

Senti

Pure Tone and Speech Audiometry,
conventional and image based in ONE instrument

Senti is modular ...

Audiometric tests

- + Pure-Tone Audiometry with bone and air conduction, talk through, external input...
- + MAGIC: Multiple-Choice Auditory Graphical Interactive Check

Speech tests in multiple languages

- + SUN: Speech understanding in Noise screening test
- + UST: Universal Speech Test adaptable to different languages
- + MATCH: Mainzer audiometric test for children - image based speech test for children from 2-5 years
- + MAUS: screening test for auditory perception and processing orders for children from 6-11 years

Additional modules available for your language and region

Product Features

- + Large and easy to use touch display of 5.0"
- + Stable desktop housing
- + USB interface to PC
- + External input modality for speech tests e.g for CD
- + Internal microphone to allow communication between patient and audiologist including the ability to monitor environmental noise and optimise measurements
- + Multiple languages support
- + Online help and information available on the device
- + Database to store 1000 patients with multiple tests per patient
- + Direct printing to label printer or via connectivity to PC
- + Updateable any time



Test Modules

- + Pure tone audiometry from 125 Hz up to 16 kHz
- + Levels from -10 to 110 dB HL in 5dB steps
- + Multiple options for accessories: HDA280, HDA200, DD45 (similar to TDH39), B71, free field, insert earphones, patient response switch
- + Bone and air conduction with masking
- + Pure tone audiometry (class 3 EN60645-1)
- + Free-Field option available

Image Based Audiometry

- + Image based audiometry and screening (MAGIC) for children > 4 years

Speech tests:

- + Multiple Speech tests (SUN, MATCH, MAUS, UST...)
- + Screening tests



New ideas in audiometry



PATH medical GmbH
Landsberger Str. 65
D-82110 Germering
Germany

Made in Germany

100366_Vero_eng_revE

Distributed By:

HEALTHWORLD INTERNATIONAL PHILIPPINES CORPORATION

0421 Quirino Ave., La Huerta, Parañaque City, Metro Manila 1700 Philippines

(+63) 2 8829-0175 | (+63) 919-065-7403 | (+63) 917-822-2938 | (+63) 998-581-2446

info@healthworldintl.com