



HEALTHWORLD
INTERNATIONAL PHILS. CORP.

World Class QUALITY Need Not Be Expensive

FDA PHILIPPINES-REGISTERED

nanoTymp®

Tympanometry simplified.

**FOR
WINDOWS**



PATH MEDICAL GmbH

Landsberger Straße 65

82110 Germering

Germany

Tel +49 89 800 76 502 / Fax +49 89 800 76 503


info@path-medical.de / www.path-medical.de

CE
0124



World Class *QUALITY* Need Not Be Expensive

PATH MEDICAL a German based company, who prides itself on developing audiological solutions, brings you the world's first PC based screening tympanometer. The nanoTymp software is compatible to most operating systems, providing tried & trusted results in a cost effective manner. *Different operating systems are compatible with the nanoTymp software - available from Q1 2020 onwards.

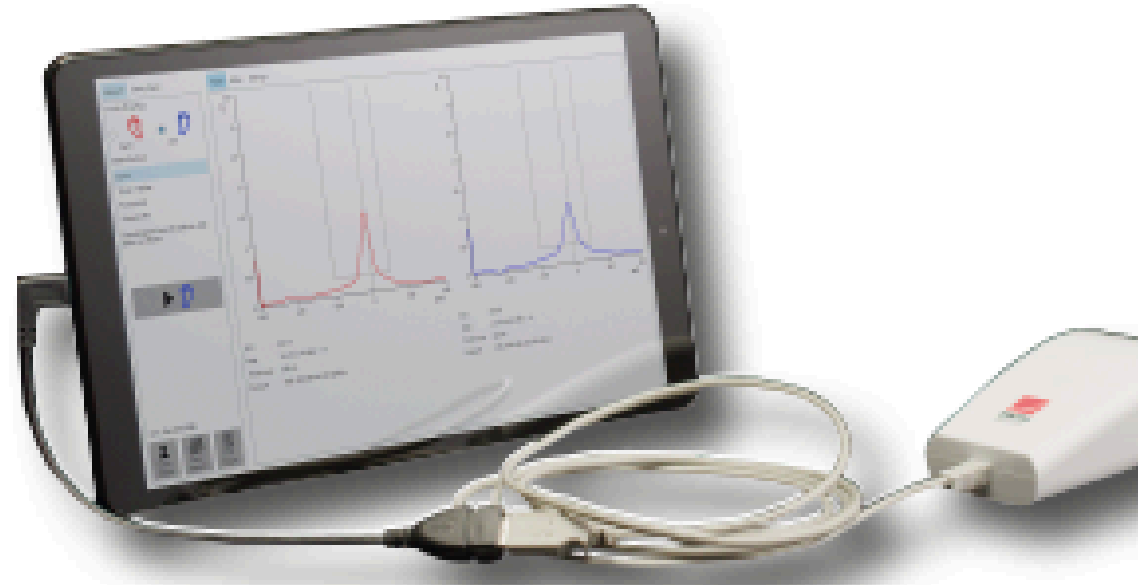

From pediatric to adult patients the nanoTymp is an invaluable tool for the verification and diagnosis of middle ear pathologies. The test procedure is easy & fast, simply fit the probe, start the test and within seconds the test results will appear on the screen. With the supplied nanoTymp software, test results & patient data can be effortlessly printed or immediately transferred via PDF to an electronic medical record.

Standard Features:

- ✓ Extra long probe cable
- ✓ 226Hz probe tone
- ✓ IPSI reflex frequencies 500 Hz, 1kHz, 2 kHz, 4 kHz, BBN
- ✓ Customized protocols
- ✓ NOAH module
- ✓ Unlimited data storage (PC dependent)


Optional features:

- 1000Hz probe tone

What's included ?

- ✓ Lanyard for hands free testing
- ✓ User manual
- ✓ Probe clip
- ✓ Calibration cavities
- ✓ nanoTymp PC software on USB stick
- ✓ Ear tip kit with probe tips

Full colour printouts can be created in seconds

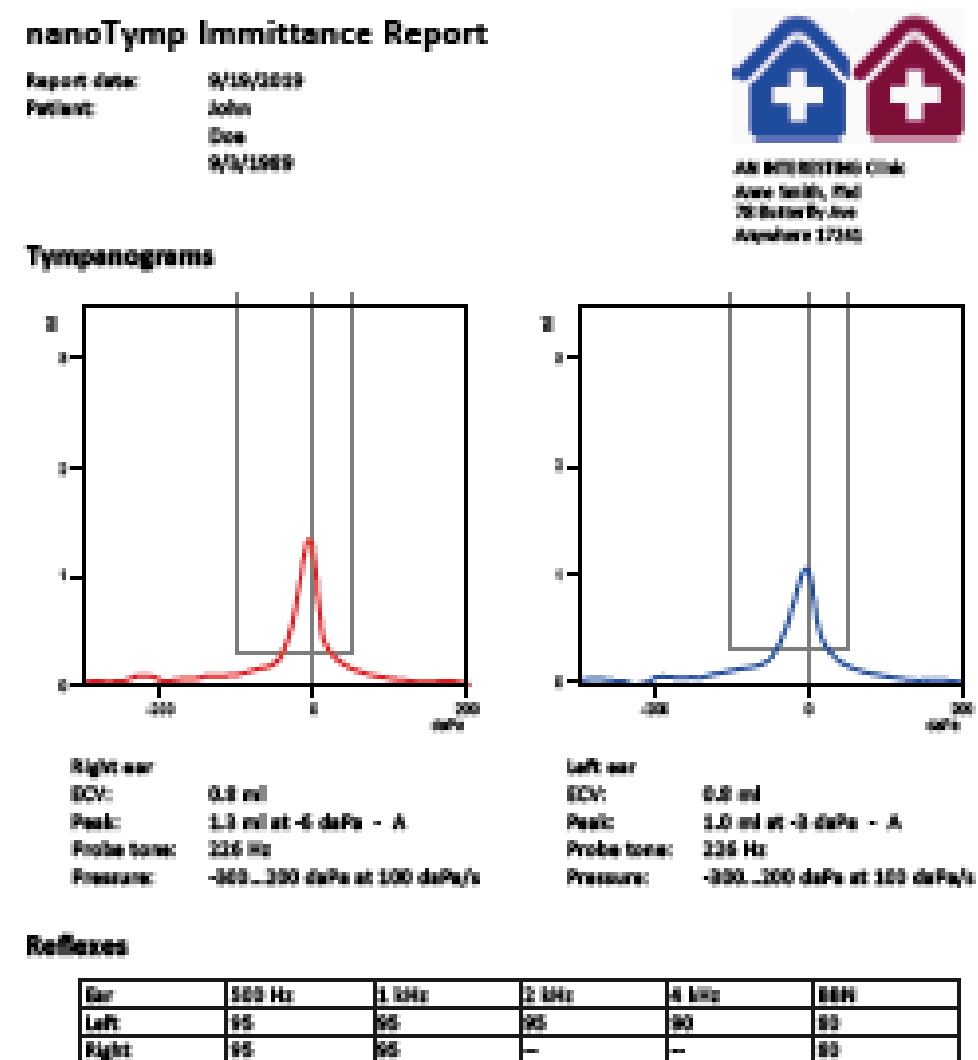

Simple USB connection: It's easy to connect nanoTymp via USB to your PC and operate it in the nanoTymp software. Unlike other Screening Tympanometers you can preview the left & right ear in real time.

Technical Specifications

Device dimensions:	150 x 72 x 55mm
Device weight:	ca. 170g
Cable length:	180 cm
Probe tone:	226 Hz ±1 % at 85.3 dB SPL ± 3 dB
Compliance range:	0 to 5 ml (±5 % or ±0, 1ml whichever is greater)
Pressure range:	- 600 to +400 daPa (±10 % or ±10 daPa whichever is greater)
Pump speed:	50, 100, 150, 200 daPa/s ± 10 daPa/s
Reflex tones (ipsi):	500 Hz, 1kHz, 2 kHz, 4 kHz, BBN
Power supply:	USB powered
Safety standards	IEC 60601-1 IEC 60601-1-2
Performance	IEC 60645-5, class 2 ANSI S3.39 type 2

Made in Germany

100937-17_Brochure_NanoTymp_EN_01 , October 2019

Distributed By:

HEALTHWORLD INTERNATIONAL PHILIPPINES CORPORATION

0421 Quirino Ave., La Huerta, Parañaque City, Metro Manila 1700 Philippines

(+63) 2 8829-0175 | (+63) 919-065-7403 | (+63) 917-822-2938 (+63) 998-581-2446

info@healthworldintl.com